



Installation Manual

Scientific Work:

Dr. Frank Götz

Editorial Work:

Dr. Frank Götz

1st Edition, Dresden, 2010

Version 2.0

Copyright © 2006 - 2010 Qualitytype AG Dresden. All rights reserved.

Production and distribution:

Qualitytype AG

Moritzburger Weg 67

D-01109 Dresden

Tel.: +49 (0)351-8838 2800

Each and every use of this product which leads to a change of the pre-defined parameters of the software is considered to be **uncertified**. Qualitytype AG does not accept any liability for improper or uncertified use of this product. Any improper or uncertified use of GenoProof® 2 is solely at the user's risk and responsibility.

GenoProof® 2 is a registered trademark of Qualitytype AG. All data in this document can be changed at any time without further notice. You are not allowed to copy or transmit this document, or parts of it, in any way for any purpose unless you have a written permission of Qualitytype AG. All product and company names are property, trademark or registered trademarks of the respective company.

Table of Contents

1.	Important Information	4
2.	System Requirements	4
3.	Setup	5
4.	Licensing.....	6
5.	Language Selection	7
6.	Getting Started.....	7

1. Important Information

Each and every use of this product which leads to a change of the pre-defined parameters of the software is considered to be **uncertified**. The Qualitytype AG does not accept any liability for improper or uncertified use of this product. Any improper or uncertified use of GenoProof® 2 is solely at the user's risk and responsibility.


GenoProof® 2 is a registered trademark of Qualitytype AG. All data in this document can be changed at any time without further notice. You are not allowed to copy or transmit this document, or parts of it, in any way for any purpose unless you have a written permission of Qualitytype AG. All product and company names are property, trademark or registered trademarks of the respective company.

2. System Requirements

- Microsoft Windows 32 Bit operating system, e.g.:
 - Microsoft Windows XP
 - Microsoft Windows Vista
 - Microsoft Windows 7
 - Microsoft Windows Server 2003
 - Microsoft Windows Server 2008
- at least 521 MB RAM
- approx. 500 MB free hard disc space (for the program without database)
- minimum screen resolution of 1024 x 768 pixels

3. Setup

1. Boot your operating system and close all active applications.
2. If there is already an older version of **GenoProof 2** installed on your computer, you must **uninstall** this version.

 **This does not apply to versions of GenoProof 1. You do not need to uninstall those versions.**

To uninstall, go to the system control (**My computer/System control**) and select **Software**. In the list, select GenoProof® 2. The uninstalling process is then started by pressing **Remove**. GenoProof® 2 is removed now. Your data including the reference data will stay unaffected. They will remain available for the next version.

 **Note:**

It might be that you are not authorized to access the system control. In this case, talk to your administrator, please.


3. Insert the CD-ROM into the optical drive and launch the desired setup program. For the desktop edition, you need to execute **GenoProof2Desktop.exe**.
4. The installation wizard will guide you. In the course of the process, the installation wizard will propose a target directory. It is possible to change the directory. However, never install more than one product in the same directory.
5. After successful installation, GenoProof® 2 can be started using the shortcut in the Windows Start menu that has been created during the installation.

4. Licensing

GenoProof® 2 can only be used, if you have a valid license key. The key can be purchased from Qualitytype AG. However, there is also the possibility to request a **free demo license key**, which enables you to test GenoProof® 2 for a limited period.

Please note that in order to generate a license key, we need to know your **system identification**. The system identification **does not contain any private information** about your system.

- When you launch the program for the first time, you will see a dialog informing you that there is **no valid license key**.

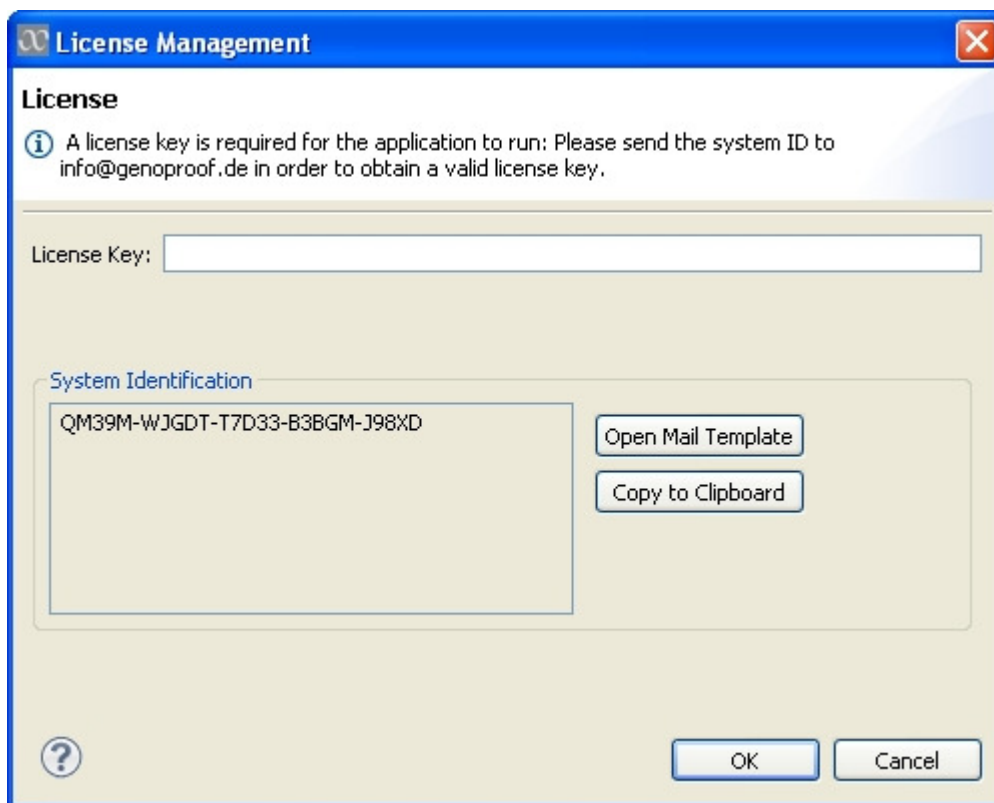

- Copy the **system identification** into an **Email** and send it to: info@genoproof.de, please.
- We will send a license key to you. Enter it into the field **License Key** and confirm it by pressing **OK**.

 **Please keep your license key safe!**

5. Language Selection

GenoProof® 2 is by default loaded with the language settings set in the **operating system** for standards and formats. However, you can determine another language that is to be used **independently** from the operating system settings.

At the moment, you can choose between the following languages:

- **German** (ISO language code: *de*)
- **English** (ISO language code: *en*)

The default language can be changed in the initialization file. To do so, proceed as follows:


1. Quit GenoProof® 2.
2. Open the **initialization file (GenoProof 2.ini)**. This file is located in the same directory like the executable file (default path: **C:\Programs\Qualitytype\GenoProof 2\Desktop**). Die initialization file has the same name like the executable file, but another file extension (**.ini**).
3. Add the following line: **-Duser.language=XX** with XX replaced by the two-place lower-case ISO-code for the referring language. The line has to be situated after the line **-Vmargs**. In addition, you can set a country format. To this end, add the line **-Duser.country=XX** with XX replaced by the two-place upper-case ISO-code for the referring country. The line has to be situated after the line **-vmargs**, too.
Save the file.



Note:

If you use a newer operating system (such as Microsoft Windows Vista), it might be that you are not authorized to access the program folder. In this case, ask your administrator to carry out the modification, please.

6. Getting Started

If you use the software for the first time, we also want to recommend you the [Getting Started Guide](#). Explanations of all functions of the program as well as the complete content of this document can be found in the program's [internal help](#). The help can be opened by pressing the  **Help** button in the toolbar of in the main window.



Qualitytype AG
Moritzburger Weg 67 | D – 01109 Dresden
Telefon +49 (0)351 – 8838 2800 | Fax +49 (0)351 – 8838 2809
www.qualitytype.de | info@genoproof.de